

THE INVACARE

Action Range



CATALOGUE



Yes, you can.®

Making Life's
Experiences Possible





Contents

Action 1R.....	4
Action 2NG	6
Action 3NG	8
Action 3NG Junior	10
Action 4NG HD	12
Action 5	14
Action Ampla.....	16



Yes, you can.®

Action 1R

A FEATURE RICH, ENTRY LEVEL CHAIR



An entry level wheelchair that doesn't compromise on functionality.

**ONLY AVAILABLE
IN NEW ZEALAND**



Action 1R

The Action 1R is an entry level wheelchair that doesn't compromise on functionality. Built using decades of in-house design experience, the Action 1R is available in transit and self-propelling versions in three sizes. The removable, swing away armrests are easy to remove and the armpad can be used in two positions.

Features



With years of experience and using tried-and-tested technology, the Action 1R is a reliable partner, the one you can count on.



The optional Eco elevating legrests provide extra comfort for users.



Optional cane holder can be mounted to securely hold a cane or crutches.



Dynamic anti-tippers provide extra safety on slopes or when backrest is reclined.

Tech Specs



Seat height
49.5cm



Seat depth
42.5cm



Total width
SW + 20cm
(self-propel)
SW + 16cm (transit)



Total length
108cm



Seat width
40.5, 45.5, 50.5cm



Max. user weight
125kg



Backrest height
43cm



Yes, you can.®

Action 2NG

SIMPLE, ROBUST AND EFFICIENT



The Action 2NG provides easy handling making it the wheelchair ideal for basic needs with quality components to give effortless rolling performance both in doors and out. A great bonus to the Action 2NG is it's cross compatibility throughout the Action range, offering familiarity and flexibility to a hugely successful Action family.

**ONLY AVAILABLE
IN AUSTRALIA**



Action 2NG

With a total weight starting from 15kg, the Invacare Action 2NG fulfills all essential functional requirements. The stable folding frame means the wheelchair is easy to handle, and being lightweight makes transportation trouble-free.

Cross-compatibility between a number of components on the Action range make this an ideal fleet product that can be easily adapted to meet individual needs. Available in self-propelled and transit versions.

Features



Decrease seat depth by 50mm for customised seat size.



The 80° footrest swings in and out for easy transfers in or out of the chair.



The Axle position can be adjusted in height by 5cm.



Lightweight foldable frame.

Tech Specs



Seat width
40.5cm, 45.5cm,
48cm



Seat depth
40-43cm



Seat height
48.5cm



Backrest height
43cm



Total width
seat width + 20cm
(self-propel)
seat width + 14.5cm
(transit)



Total length
108cm



Max. user weight
125kg



Transport weight
12.1kg (transit)
8.9kg (self-propel)



Total weight unloaded
14.9kg (transit)
16.2kg (self-propel)



Turning radius
87cm



Max. safe slope
9°



Crash Test Approved
Fulfills the requirements of ISO 7176-19



Yes, you can.®

Action 3NG

GREATER DRIVING PERFORMANCE



The Action 3NG manual wheelchair offers a sleek and stylish chair for everyday use. Delivering outstanding drive performance with new lightweight options for individuals with a more active lifestyle. Highly functional yet elegantly designed the Action 3NG wheelchair is a lightweight chair offering a vast range of options.

Action 3NG

Adaptable and modular, the Action®3NG accommodates progressive changes in user needs. It is also simple to refurbish, with a wide range of options.

With a low transfer weight the Action®3NG is easy to fold and to dismantle, making it easily transportable in a small car.



Features



Brake and Brake extension. Fold the brakehandle for easier side transfers.



Action 3NG Lever Drive.



20 rear wheel positions.



Magnesium castor attachment.

Tech Specs



Seat width
40.5cm, 45cm,
50.5cm



Seat depth
40.5cm, 45cm



Seat height
48.5cm



Backrest height
43cm



Total width
seat width + 16cm
(transit)
seat width + 19cm
(self-propel)



Total length
108cm



Max. user weight
125kg



Transport weight
10.3kg (transit)
7.1kg (self-propel)



Total weight unloaded
13.5kg (transit)
13kg (self-propel)



Turning radius
87cm



Max. safe slope
9°



Yes, you can.®

Action[®] 3 Junior

YOUR JOURNEY, YOUR WAY



Invacare's Action 3 Junior is a lightweight, foldable paediatric wheelchair designed for children aged between 3 and 15 years. Action 3 Junior has been developed to match the individual needs of the child and can grow as they grow. The Action 3 Junior comes in attractive, modern styles, designed to appeal to kids.



Action 3 Junior

As children grow and develop, a key feature of a children's wheelchair is that it can be adjusted to accommodate the child's individual growth and needs. With the Invacare Action 3 Junior, the width of the seat and backrest can be adjusted, as can the height of the arm and footrests.

A wide range of options are available on the Action 3 Junior to meet every child's needs in terms of comfort and positioning.

Features and Options



Innovative telescopic cross brace
allows the seat to be adjusted in width and depth. With the seat and backrest upholstery also allowing for quick and easy width adjustability, there is no need for additional parts.



Upholstery adjustment
As the chassis is made wider, the upholstery on the back and on the seat of the chair can be made wider with Velcro flaps.



Front stabilizer kit
A front stabilizer kit is available to increase forward stability for active children who are likely to climb all over the wheelchair and be especially energetic when propelling the wheelchair.



Adjustable seat height and balance
For improved driving performance and different height options, the rear wheel can be mounted in 20 different positions.

Tech Specs



Seat width
20.5 - 38cm



Seat depth
20 - 37.5cm



Seat angle
1.5° / 5° / 8.5°



Seat height
41 / 43.5 / 46cm



Backrest angle
0° / 7° / 0° - 30°



Backrest height
35 - 43cm



Front frame angle
80° / 90° /
0° - 80°



Length of footrests
12.5 - 45cm



Armrest height
15-25 / 18-28cm



Driving unit width
Seat width +
19/17cm



Total folded width
31.5/27*cm



Total length including footrests
74 - 100cm



Total length excluding footrests
51 - 77cm



Total height
87 - 113cm



Max. safe slope
12°



Turning radius
75cm



Max. user weight
80 kg



Total weight from
15 / 14.5' kg



Mass of the heaviest part, from
8.5 / 12' kg



1. Transit version



Yes, you can.®

Action 4NG HD

STABILITY MEETS STYLE



The Action 4NG HD has outstanding functionality with a user weight of 160kg, offering flexibility in a chair to suit almost everyone.

Action 4NG HD

A comfortable, foldable, lightweight wheelchair, the Action 4NG is ideally suited for everyday use.

A key feature of this next-generation model is the innovative Invacare hallmark of cross-compatibility; a shared platform with other models in the Action NG family means that users and advisors can easily update or refurbish their wheelchairs and enjoy maximum opportunities for personalised adjustment.



Features



Tech Specs



Seat height
51cm



Max. safe slope
9°



Total weight unloaded
14kg



Seat width
50.5cm, 55.5cm,
60.5cm



Turning radius
85cm



Transport weight
8kg



Seat depth
45-50cm



Total width
seat width + 20cm



Max. user weight
160kg



Backrest height
43cm



Total length
117cm



Crash Test Approved
Fulfills the requirements of ISO 7176-19



Yes, you can.®

Action 5

THE PERFECT FIT FOR YOUR LIFESTYLE



The Invacare Action⁵ has been designed for the more active user looking for an affordable chair with good drive performance and the ability to fold easily for handling, storage and transportation purposes. With its two patents pending and lightweight construction, the new Action⁵ gives the user a chair they can rely on.

**ONLY AVAILABLE
IN AUSTRALIA**

Action 5



Designed specifically for users with an active lifestyle, the Invacare Action 5 is a highly adjustable, stable chair with a unique and innovative folding mechanism that allows the backrest to fit neatly on top of the seat.


The Action 5 features a lightweight frame with step-less adjustability, meaning very fine adjustments can be made. In addition many of the chair's features can be adjusted while the user remains seated. Swing away and removable armrests allow the Action 5 to be reconfigured for a variety of different situations.

The Action 5 is available in a wide range of attractive colours as well as three accent colours for smaller parts like the castors and rear wheel attachments.

Features



Tension and depth adjustments can be made via easy-to-use straps.

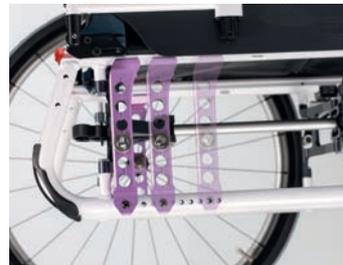

100mm centre of gravity adjustments can be made in 25mm increments which can even be fine-tuned to 12mm.



Standard adjustments of 0-2° available with the option of 3 & 4° adjustments if required.



A wide range of options and accessories are available to give your Action 5 a personal touch.

Tech Specs



Seat width

30.5 / 33 / 35.5 /
38 / 40.5 / 43 /
45.5 / 48 / 50.5 /
53 cm



Seat depth

35 / 37.5 / 40 /
42.5 / 45 / 47.5 /
50 cm



Seat height

front: 36 - 53.5cm
rear: 36 - 50cm
(increments of
2.5cm)



Total width

SW +16 cm To
SW+ 25 cm



Total width folded

32 -38 cm



Total length

84.5 -120.5 cm



Max. user weight
100 & 130 kg



Transport weight
w/out rear wheels
8.3 kg



Total weight
unloaded
from 12.3kg



Turning radius
105 -138 cm



Max. safe slope
7°



Crash Test
Approved
Fulfills the
requirements
of ISO 7176-19



Yes, you can.®

Action[®] Ampla[™]

PROVIDING INDEPENDENCE FOR AN ACTIVE LIFESTYLE



TAKING BARIATRIC WHEELCHAIRS TO THE NEXT LEVEL

INNOVATIVE | AGILE | AWARD WINNING



reddot design award
winner 2019





Taking bariatric wheelchairs to the next level

The Invacare **Action® Ampla™** has been designed and developed specifically for Plus-Size users, taking into consideration body shapes, varying activity levels and weight limit requirements. But it doesn't stop there! The sleek and minimalist wheelchair, with its modern appearance, provides independence and support for a comfortable and active lifestyle.



Award winning

Already an award-winning chair, the agile and stylish Invacare **Action® Ampla™** has won a Red Dot Award for Product Design. This award is not only the largest product competition in the world, but also one of the most prestigious design awards. The **Action® Ampla™** has also won the Italian Platinum A'Design Award and the Gold German Design Award.



Yes, you can.®

5

KEY FEATURES

1

Inbuilt transportation solution

The transportation bearings on the **Action® Ampla™** have been cleverly designed to take the weight of the wheelchair, allowing it to easily go through doorways or narrow spaces:

1. Fold down the backrest
2. Lift the legrest section of the wheelchair
3. Tilt wheelchair onto the transportation bearings

Whilst in this position, the wheelchair can be moved easily.

1.



2.



3.



5

Multiple active rear wheel positions

- ▶ Very easy to propel
- ▶ Users can go further, using less energy



2

Ergonomic backrest

- ▶ Maximum access to handrims
- ▶ No obstructions
- ▶ More freedom of movement
- ▶ Superb driving performance and manoeuvrability



3

Height adjustable transfer levers

- ▶ Support for entering and leaving the chair
- ▶ Can be turned inwards or outwards
- ▶ Can be lowered to allow table access



4

Personalised accent

- ▶ Three stylish colour choices





Yes, you can.®

Innovative Backrest Design



The design of the wheelchair allows for a perfect fit for both Plus-Size male users, who tend to have more of an "apple" body-shape and Plus-Size female users, who tend to have a "pear" body-shape. This is thanks to the innovative backrest on the **Action® Ampla™** that has been designed to accommodate the anatomy of both body-shapes whilst offering support and comfort to the user.

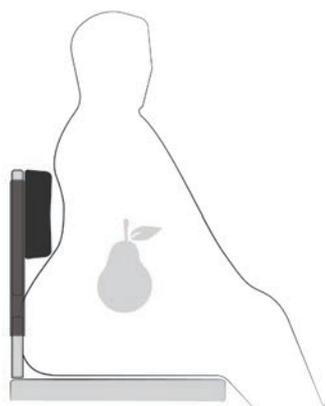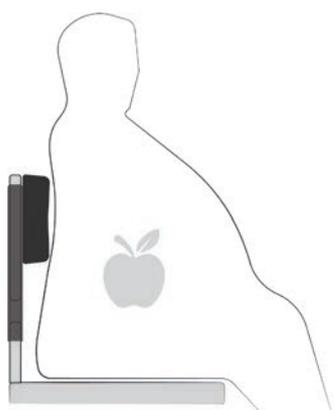

Main features

- ▶ Inbuilt transportation solution
- ▶ Ergonomic backrest
- ▶ Height adjustable transfer lever
- ▶ Personalised accent
- ▶ Multiple active rear wheel positions

Additional features


- ▶ An elastic, breathable privacy cover hides the gluteal shelf
- ▶ Push to lock brakes are nicely integrated into the side guard for accessibility

Tech Specs


For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

 <p>Seat width 50.5cm, 55.5 cm, 60.5cm, 65.5cm, 70.5cm, 75.5cm</p>	 <p>Seat depth 46-51cm</p>	 <p>Seat height 38 - 45.5cm in steps of 2.5cm</p>	 <p>Backrest height 45 - 52.5cm</p>	 <p>Total width seat width + 20cm (self-propel) seat width + 14.5cm (transit)</p>	 <p>Total length 108cm</p>
 <p>Max. user weight 300kg</p>	 <p>Transport weight From 19.5kg</p>	 <p>Total weight From 27kg</p>	 <p>Turning radius 142 - 166cm</p>	 <p>Max. safe slope 7°</p>	 <p>Crash Test Approved Fulfills the requirements of ISO 7176-19</p>

Frame Colours

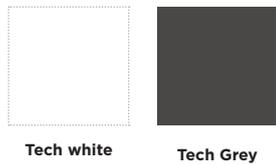

Accent Colours

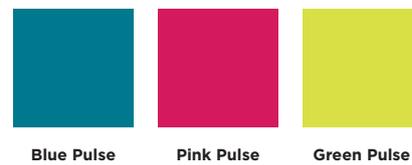

Please note: colours may vary slightly from those displayed here.



FOR FURTHER INFORMATION:

New Zealand: Freephone 0800 468 222, email sales@invacare.co.nz or visit www.invacare.co.nz

Australia: Freephone 1800 460 460, email orders@invacare.com.au or visit www.invacare.com.au

All rights reserved. All information quoted is believed to be correct at time of publication.
Invacare reserves the right to alter product specifications without prior consent. **05/2024**



Yes, you can.®