

Azalea®

MAKES LIFE COMFORTABLE



OFFERS EXCELLENT SEATING COMFORT AND
MANY ADJUSTMENT POSSIBILITIES



Yes, you can.®



The **Rea Azalea** is the natural choice for the moderately active client as well as for the more passive one. The comfortable seating system allows the client to alter the seat position, ensuring maximum comfort at all times, thus combined with excellent driving characteristics makes **Azalea** the ideal choice. The **Azalea** is available as a self-propelling wheelchair with large rear wheels or an alternative transit version that keeps the same frame but with smaller rear wheels. The Flex 3 backrest is multi-adjustable in height, width and shape for different body forms.

Tilt & recline makes life comfortable

The **Rea Azalea** offers 25° seat tilt and 30° backrest recline. Both tilt and recline can be easily adjusted manually or electronically by the assistant and also by the client with the electronic option. The backrest has a biomechanical positioning hinge, which follows the user's natural back movement when adjusting the backrest angle.

Combined with the seat-tilt function, the DSS (Dual Stability System) helps to prevent Decubitus thanks to its non-friction body repositioning.



Key features and benefits



▲ Flex 3, backrest plate

- The backrest plate is a strong and sturdy plate offering the client stability whether upright or in tilt/recline.
- This backrest plate offers height adjustment for the taller clients and width adjustment for wider clients (wider or taller back cushions may be required to accommodate the change in size).
- Standard feature on stock models in Australia and New Zealand.



▲ Flex 3, tension adjustable backrest

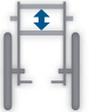
- The multi adjustable backrest can be adjusted in height and in width for different body shapes.
- You can even adjust it so it is wider at the top and narrower at the bottom or vice versa, for those clients that have wider shoulders or hips.



▲ Other adaptations and accessories

- Adjustable backrest height, seat width and seat depth, the chair can be easily adjusted to the persons exact specification.
- Individual width adjustments for the arm rests and leg rests for individual client needs.
- Same chassis for all rear wheel size, reducing the number of parts required to accommodate different needs.
- The client or carer can easily change the seating position with the power tilt and recline option.

Tech Specs

 <p>Seat width 39-44/ 44-49/ 49-54cm</p>	 <p>Seat depth 43 - 50cm</p>	 <p>Seat height 40/45cm</p>	 <p>Backrest height Backrest plate 56-79cm Tension adjustable backrest: 54-71.5cm</p>	 <p>Armrest height 24-36cm</p>	 <p>Footrest length 33 - 62.7cm</p>	 <p>Seat Angle 1 -25°</p>	 <p>Backrest Angle 0-30°</p>
 <p>Total width seat width + 22cm</p>	 <p>Total height 90-125cm</p>	 <p>Total length 102cm (inc footrests)</p>	 <p>Total weight 34kg (SW 440 mm)</p>	 <p>Max. user weight 135kg</p>	 <p>Transport weight 16.5kg</p>	 <p>Crash Tesst Approved ISO 7176-19</p>	 <p>Max. safe slope 7°</p>

Frame Colour



Pearl Grey



Azurite Blue



Sand



Happy Red



Electric Green



Yes, you can.®



FOR FURTHER INFORMATION:

New Zealand: Freephone 0800 468 222, email sales@invacare.co.nz or visit www.invacare.co.nz

Australia: Freephone 1800 460 460, email orders@invacare.com.au or visit www.invacare.com.au

Asia: Phone +61 2 8839 5330, email anzexportsales@invacare.com or visit www.invacare.com

Thailand: Phone +66 0 2821 5515 or visit www.invacare.com



Yes, you can.[®]